Quantum Blue[®] Calprotectin

Quantitative Lateral Flow Assay

LF-CAL25

25 tests

Revision date: 2015-02-09
Version 2

ENGLISH

INTENDED USE

BÜHLMANN Quantum Blue® Calprotectin assay is an immunoassay designed for the quantitative determination of Calprotectin concentrations in human stool samples (ref. 1-6) in combination with the BÜHLMANN Quantum Blue® Reader.

For professional use only*.

*Canada, Taiwan: For laboratory use only.

PRINCIPLE OF THE ASSAY

The test is designed for the selective measurement of Calprotectin antigen by sandwich immunoassay. A monoclonal capture antibody (mAb) being highly specific for Calprotectin is coated onto the test membrane. A second monoclonal detection antibody conjugated to gold colloids is deposited onto the conjugate release pad and released into the reaction system after addition of the extracted and diluted stool sample. The Calprotectin/anti-Calprotectin gold conjugate binds to the anti-Calprotectin antibody coated on the test membrane (test line; test band) and the remaining free anti-Calprotectin gold conjugate binds to the goat anti-mouse antibody coated on the test membrane (control line; control band). The signal intensities of the test line and the control line are measured quantitatively by the BÜHLMANN Quantum Blue® Reader.

REAGENTS SUPPLIED AND PREPARATION

Reagents	Quantity	Code	Comments
Test Cassette	25 cartridges	B-CAL-TC	vacuum-sealed in a foil bag
Extraction Buffer	1 bottle 125 mL	B-CAL-EX	Ready to use
RFID Chip Card	1	B-LFCAL-RCC	White plastic card

Table 1

REAGENTS & MATERIAL SUPPLIED SUPPLEMENTARY

Fecal Extraction Devices

Fecal Extraction Devices described below are not delivered with the kit and either of them **has to be ordered** with the kit.

CALEX® Cap Device	50 or 200 tubes filled with 5 mL extraction buffer, ready to use	B-CALEX-C50 B-CALEX-C200
Smart-Prep	50 tubes consisting of spatulas and base caps	B-CAL-RD
ScheBo® Quick-Prep™	50 tubes consisting of tube, cone & dosing tip. Filled with 1.3 mL Extraction Buffer, ready to use	B-CAL-SO50

Table 2

Controls

Controls for BÜHLMANN Quantum Blue® LF-CAL25 can be ordered separately.

Controls Low* / High*	2 vials, 0.5 mL	B-CAL-CON	Ready to use
-----------------------	--------------------	-----------	--------------

Table 3

* Lot specific concentration of Calprotectin

STORAGE AND SHELF LIFE OF REAGENTS

All kit components are stable at 2-8°C until the expiration date printed on the labels.

PRECAUTIONS

SAFETY PRECAUTIONS

- None of the reagents of this test contains components of human origin.
- Patient specimens should be handled as if capable of transmitting infections and should be handled in accordance with Good Laboratory Practice (GLP) using appropriate precautions.
- Unused solution should be disposed of according to local State and Federal regulations.

TECHNICAL PRECAUTIONS

Kit components

- All reagents and test samples must be equilibrated to room temperature (18-28°C) before starting the assay.
- Components must not be used after the expiry date printed on the labels.
- Do not mix different lots of reagents.
- Test Cassettes cannot be re-used.

Test Procedure

- Read carefully the instructions prior to carrying out the test. Test performance will be adversely affected, if reagents are incorrectly diluted, handled or stored under conditions other than those as detailed in this instruction for use.
- The Quantum Blue® Reader must be switched on and programmed for the Quantum Blue® Calprotectin assay (CAL_0 or CAL_720) before starting the assay (see Quantum Blue® Reader Manual).
- Use the white RFID Chip Card in order to change lot-specific test parameters.
- Patient samples that are not properly handled may cause inaccurate results.
- In order to receive reliable and quantitative results it is important to homogenize the stool sample entirely in the Extraction Device.
- Diluted samples should be used within several hours and cannot be stored for a longer time period.
- With BÜHLMANN Smart Prep and ScheBo® Quick-Prep™, it is important to centrifuge the extracts before storage. Centrifuge the tubes for 5 minutes at 3'000 x g. After centrifugation the supernatant must be transferred into a fresh storage tube. With CALEX® Cap Device you can store the extracts directly with the device. For re-use/re-measurement of the extracts see step 2 under the chapter Assay Procedure.

MATERIALS REQUIRED BUT NOT PROVIDED

- Vortex mixer for stool extraction
- Precision pipettes with disposable tips: 10-200 µl and 1 ml
- Centrifuge
- 5 ml polypropylene or polystyrene tubes for dilution of the extracts
- Timer (optional)
- Quantum Blue® Reader available from BÜHLMANN (order code: BI-POCTR-ABS)
- Soft tissues or blotting paper

SPECIMEN COLLECTION AND STORAGE

If the Extraction Devices are used (refer to Table 2), less than 1 g of native stool will be required.

Collect stool samples into plain tubes. The samples can be stored refrigerated at 2-8°C for at least 6 days.

The extracts are stable for at least 7 days at 2-8°C and for at least 24 months at ≤ -20°C.

Important: The sample must be collected without any chemical or biological additives in the collection device.

ASSAY PROCEDURE

The Assay Procedure consists of three steps:

1. Extraction of stool samples:

The extraction is described in the instruction for use delivered with the respective extraction devices.

2. Sample processing:

- Smart Prep or ScheBo® Quick Prep™: After extraction, let the stool extract settle for 10 minutes. Dilute the supernatant 1: 16 with extraction buffer (e.g. 20 µL extract and 300 µL extraction buffer) and mix well. Let the samples equilibrate for at least 5 minutes at 18-28°C prior to proceeding to the next step (step no. 3).

- CALEX® Cap Device: After extraction, let the stool extract settle for 10 minutes with the white head of the device down. Unscrew the blue cap and dilute the supernatant 1:1.6 with extraction buffer (e.g. 200 µL extract and 120 µL extraction buffer) and mix well. Let the samples equilibrate for at least 5 minutes at 18-28°C prior to proceeding to the next step (step no. 3).

3. Lateral Flow Assay Procedure and Readout:

There are two alternative methods available on the Quantum Blue® Reader: <CAL_720> and <CAL_0>. Select one of these methods before starting the experiments.

Load the lot specific parameters from the red RFID Chip Card.

3.1. Method <CAL_720> with internal timer

- Load the test cassette onto the test cassette holder of the Reader.
- Add 60 µl of diluted stool extract onto the sample loading port of the test cassette.
- Close the cassette holder and start the measurement by pressing the start button.
- The scan starts automatically after 12 minutes (720 seconds).

3.2. Method <CAL_0> without internal timer

- Add 60 µl of diluted stool extract onto the sample loading port of the test cassette.
- Incubate for 12 minutes +/- 1 minute (set a timer manually).
- Load the test cassette onto the test cassette holder of the Reader.
- Scan the cassette with the Quantum Blue® Reader by pressing the start (<ENTER>) button immediately.

Remark: Please refer to the Quantum Blue® Reader Manual to learn about the basic functions and how to initialize and operate the Reader, especially how to select test methods, and how to load lot specific parameters from the RFID Chip Card in order to get the samples measured.

QUALITY CONTROL

- If the precision of the assay does not correlate with the established limits and repetition excludes errors in technique, check the following issues: i) pipetting, temperature controlling and timing devices ii) expiration

dates of reagents and iii) storage and incubation conditions.

- Result of the self-test of the Quantum Blue® Reader performed at startup of the instrument has to be valid.
-

VALIDATION OF RESULTS

- For a valid test result, the Control Line (C) must be visible in any case (see Figures 1A and 1B). It is used as functional test control only and cannot be used for the interpretation of the Test Line (T). If the Test Line (T) is not detectable after 12 minutes of incubation time (Figure 1A), the concentration of Calprotectin present in the stool sample is below the detection limit. If a Test Line (T) is detectable after 12 minutes of incubation time (Figure 1B), the Calprotectin concentration present in the stool sample is calculated by the Quantum Blue® Reader.
 - If only the Test Line (T) is detectable after 12 minutes of incubation time (Figure 1C), the test result is invalid and the Calprotectin assay has to be repeated using another test Cassette.
 - If neither the Control Line (C) nor the Test Line (T) are detectable after 12 minutes of incubation time (Figure 1D), the test result is invalid and the Calprotectin assay has to be repeated using another test Cassette.
 - As the Quantum Blue® Reader allows a quantitative evaluation of the Test (T) and Control (C) Lines, an additional validity check of the Control Line (C) is undertaken. If the signal intensity of the Control Line (C) is below a threshold after 12 minutes of incubation time, the test result is also invalid and the Calprotectin High Range assay has to be repeated using another test Cassette.
-

STANDARDIZATION AND INTERPRETATION OF RESULTS

- The Lateral Flow Assay is calibrated with the BÜHLMANN fCAL™ ELISA (order code: EK-CAL).
 - The BÜHLMANN Quantum Blue® Reader uses a lot-specific standard curve to calculate the Calprotectin concentration. This lot-specific standard curve is generated with the median values (n≥10 measurements each) from at least 10 calibration points obtained from different stool samples with known Calprotectin concentrations. The assay range is between 30 and 300 µg/g.
 - For quantitative measurements, unknown samples reading above 300 µg/g can be re-tested in the BÜHLMANN Quantum Blue® Calprotectin assay (order code: LF-CHR25).
-

LIMITATIONS

- The reagents supplied with this kit are optimized to measure human Calprotectin in extracted stool samples.
- Fecal Calprotectin values can be used as supplementary data available to the physician in establishing a diagnosis.
- A Calprotectin test result obtained with the BÜHLMANN Quantum Blue® Calprotectin rapid assay can be confirmed by measuring other biomarkers (such as C-Reactive Protein) or by repeating the Calprotectin test (with BÜHLMANN fCAL™ ELISA or using a second test cassette) and comparing to the CDAI score and the clinical status of the patient, before any clinical and/or therapeutic decision is made.
- During therapy follow-up, a Calprotectin test result obtained with the BÜHLMANN Quantum Blue® Calprotectin rapid assay can be interpreted by

comparing the result with previous ones and in the context of the patient's clinical picture.

EXPECTED VALUES

Estimation of fecal Calprotectin is a reliable and simple way to distinguish organic from functional gastrointestinal diseases.

The following data were established with the BÜHLMANN fCAL™ (order code: EK-CAL). In a clinical study, 401 symptomatic consecutive patients scheduled for colonoscopy were investigated (Table 19).

Endoscopy examination showed 273 patients with functional diseases, whereas 128 patients had various organic diseases (colitis, Crohn's, ulceration, diverticulitis, polyps, adenomas, cancer, or infectious diseases).

A clinical cut-off at 50 µg/g was determined by ROC curve analysis (AUC: 0.935). By applying this cut-off for differentiation between organic and functional diseases, a clinical sensitivity of 84.4% and a specificity of 94.5% can be reached. The suggested cut-off level for adults (50 µg/g) can be used for children aged from 4 to 17 years regardless of sex (ref.7). Although fecal Calprotectin levels from adults and children are reported comparable in general, levels of newborns can be significantly increased (ref.6).

These data support the following recommendation for interpretation of results:

Values below 50 µg/g:

Calprotectin values <50 µg/g are not indicative of inflammation in the gastrointestinal tract. Patient-samples with low levels are likely not to be in need of further invasive procedures to determine the inflammation cause (ref. 4).

Elevated values between 50 and 200 µg/g:

Calprotectin values between 50 and 200 µg/g can represent mild organic disease such as inflammation caused by NSAIDs, mild diverticulitis and IBD in remission phase. The low inflammatory response shown within this range may suggest repeating the measurement and performing further investigations.

Elevated values above 200 µg/g:

Calprotectin values >200 µg/g are indicative of active organic disease with inflammation in the gastrointestinal tract. Appropriate further investigative and curative procedures by specialists are suggested.

PERFORMANCE CHARACTERISTICS

Limit of Blank (LoB): <7 µg/g calprotectin. The LoB has been established in accordance with CLSI protocol EP17-A in three independent runs using three different lots of Test Cartridges with 60 blank values in total by using Extraction Buffer as a sample.

Limit of Detection (LoD): 15 µg/g calprotectin. The LoD has been established with two stool extracts with concentrations of 14.5 and 28.5 µg/g. The samples were measured with three different lots of Test Cartridges in three independent runs of 10 replicates each, set up within 10 minutes. The averaged SD values were determined and the LoD has therefrom been calculated in accordance with CLSI protocol EP17-A.

Limit of Quantification (LoQ):

Lower LoQ: ≤30 µg/g calprotectin

Upper LoQ: ≥300 µg/g calprotectin

The LoQ has been established with seven stool extracts at concentrations between 14.5 and 622 µg/g. The samples were measured with three different lots of Test Cassettes in three independent runs of 10-20 replicates each, set up within 10 minutes and averaged. The resulting precision profile is shown in Figure 2. The limit of quantification corresponds to the concentration of Calprotectin with an imprecision below 25% CV allowing a quantitative measurement within the range from 30 (lower LoQ) to 300 µg/g (upper LoQ).

Linearity: Four stool samples with elevated Calprotectin concentrations were extracted according to the assay procedure. The extracts were diluted with Extraction Buffer and each dilution was subsequently assayed on two test Cassettes according to the assay procedure and averaged for each test point. The results showed linearity within the indicated measuring range of 30 to 300 µg/g of the Quantum Blue® Calprotectin assay for all 4 samples ($R^2 = 0.986 - 0.995$). One example is shown in (Figure 3). The average difference between measured and expected concentrations was **2%** corresponding to a recovery of 98%. The recovery varied between 84-114 %.

Recovery: 92%. Four stool extracts containing low calprotectin concentrations were spiked with different amounts from a highly concentrated calprotectin stock solution. Each sample was measured before and after spiking on three test cartridges according to the assay procedure. The three replicates were averaged and the mean values are presented in Table 21. The recovery varied 83-105 %.

Repeatability: 19% CV. The repeatability of the Quantum Blue® Calprotectin assay was calculated from 2 extracted stool samples containing 49 and 138 µg/g calprotectin. Each sample was tested according to the assay procedure in one run of 20 replicates set up within 10 minutes. The mean values of three different lots of Test Cassettes are presented in (Table 20). The repeatability varied between 15.5-22.2 % CV.

Inter-lot Precision: 22% CV. The inter-lot precision of the Quantum Blue® Calprotectin assay was calculated from 4 extracted stool samples containing between 30 and 244 µg/g calprotectin. Each sample was tested according to the assay procedure in three independent runs of 10-20 replicates each, set up within 10 minutes using three different lots of Test Cassettes. The inter-lot data of the 4 samples varied from 20.2 to 24.5% CV.

Method Comparison: $R^2 = 0.806$; $y=1.02x-11.1$ µg/g 101 stool extracts from subjects undergoing IBD vs. IBS screening exhibiting calprotectin concentrations within the indicated measuring range of the Quantum Blue® Calprotectin assay were analyzed according to the assay procedure using a data set originated from Test Cassettes of six different lots and compared with the values obtained by the BÜHLMANN fCAL™ (order code: EK-CAL). The correlation data are illustrated in Figure 4.

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Table 20 Precision (Repeatability)

Lot	Sample	Mean [µg/g]	SD [µg/g]	CV [%]
A	S1	47.8	10.6	22.2
A	S2	141.3	22.0	15.5
B	S1	50.1	10.0	20.1
B	S2	145.5	23.5	16.2
C	S1	50.0	10.3	20.6
C	S2	146.6	31.5	21.7
Mean				19.4

Table 21 Spiking Recovery

Sample	µg/g	spiked with [µg/g]	Observed [µg/g]	Expected [µg/g]	O/E [%]
P1	53	25	62	78	79
		50	154	103	150
		100	143	153	94
		200	239	253	95
Mean					105
P2	81	25	107	106	101
		50	120	131	91
		100	162	181	89
		200	265	281	94
Mean					94
P3	56	25	77	81	94
		50	120	106	113
		100	110	156	71
		200	186	256	73
Mean					88
P4	160	25	153	185	82
		50	191	210	91
		100	190	260	73
		200	305	360	85
Mean					83
Mean					92

Figure 2

Precision Profile

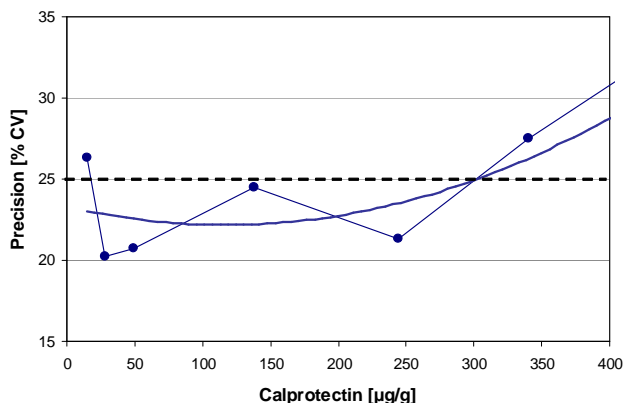
Figure 3:

Linearity Plot


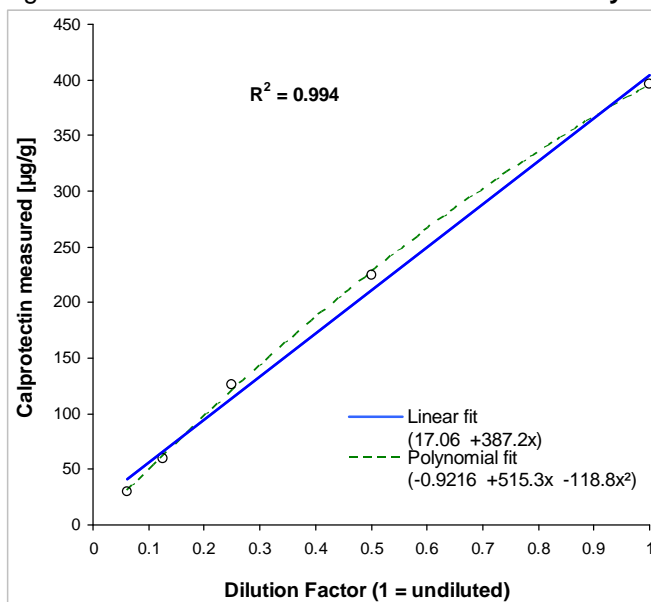
Figure 4:

Method Comparison

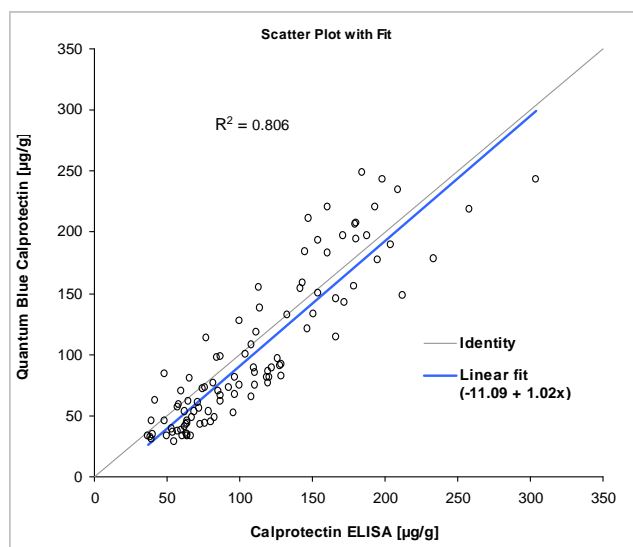REFERENCES/ LITERATURREFERENZEN/ REFERENCES/ RIFERIMENTI/ REFERENCIAS

1. Tibble JA et al.: *Use of surrogate markers of inflammation and Rome criteria to distinguish organic from nonorganic intestinal disease.* Gastroenterol 123, 450-460 (2002)
2. Fagerhol MK: *Calprotectin, a faecal marker of organic gastrointestinal abnormality.* Lancet 356, 1783-4 (2000)
3. Hessian PA and Fisher L: *The heterodimeric complex of MRP-8 (S100A8) and MRP-14 (S100A9). Antibody recognition, epitope definition and the implications for structure.* Eur J Biochem 268, 353-63 (2001).
4. Tschanguizi et al.: *Stool Calprotectin as a Marker of inflammation.* Poster at UEGW, 2007
5. Tibble JA et al.: *A simple method for assessing intestinal inflammation in Crohn's disease.* Gut 47, 506-513 (2000).
6. Konikoff MR and Denson LA: *Role of fecal calprotectin as a biomarker of intestinal inflammation in inflammatory bowel disease.* Inflamm Bowel Dis 12(6), 524-34 (2006)
7. Fagerberg UL et al.: *Colorectal inflammation is well predicted by fecal calprotectin in children with gastrointestinal symptoms.* J Pediatr Gastroenterol Nutr 40, 450-5 (2005)
8. Jahnsen J, Røseth AG, Aadland E. *Measurement of calprotectin in faeces.* Tidsskr Nor Lægeforen 128, 743-5 (2008)

SYMBOLS/ SYMBOLE/ SYMBOLES/ SIMBOLI/ SIMBOLOS

Symbol	Explanation
	Use By Verwendbar bis Utiliser jusqu'au Utilizzare entro Fecha de caducidad
	Catalogue number Bestellnummer Référence du catalogue Numero di catalogo Número de catálogo
	Batch code Chargenbezeichnung Code du lot Codice del lotto Codigo de lote
	Contains sufficient for <n> tests Ausreichend für "n" Ansätze Contenu suffisant pour „n“ tests Contenuto sufficiente per „n“ saggi Contenido suficiente para <n> ensayos
	Consult Instructions for Use- Gebrauchsanweisung beachten Consulter le mode d'emploi Consultare le istruzioni per l'uso Consulte las instrucciones de uso

Symbol	Explanation
	Temperature limitation Zulässiger Temperaturbereich Limites de température Limiti di temperatura Limite de temperatura
	Test Cassette Testkassette Cassette de testi Cartuccia di rilevazione Cartucho de prueba
	Extraction Buffer Extraktionspuffer Tampon d'extraction Tampone di estrazione Tampón de extracción
	<i>In Vitro</i> Diagnostic Medical Device <i>In Vitro</i> Diagnostikum Dispositif médical de diagnostic <i>in vitro</i> Dispositivo medico-diagnostico <i>in vitro</i> Producto sanitario para diagnóstico <i>in vitro</i>
RCC	RFID Chip card RFID Chipkarte RFID Carte à puce RFID Carta chip RFID Tarjeta chip

Printing Date
2015-02-10

Not for sale in the USA, due to current patent restrictions