

Quantum Blue® Adalimumab

trough level rapid test

adapt therapy
without delay

- First quantitative rapid test for Adalimumab monitoring
- Time to result 15 minutes
- Measuring Range 1.3-35 µg/mL
- High agreement to reference ELISAs

Health Canada Licence: 101776


Quantum Blue® Adalimumab rapid test

Key features

- The first quantitative rapid test for adalimumab
- Measuring range 1.3 -35 µg/mL
- Covers the range of the relevant 5 to 12 µg/mL for trough levels*
- High correlation and agreement to the reference ELISA assays

* Mitrev et al. 2017

Correlations to ELISA assay:

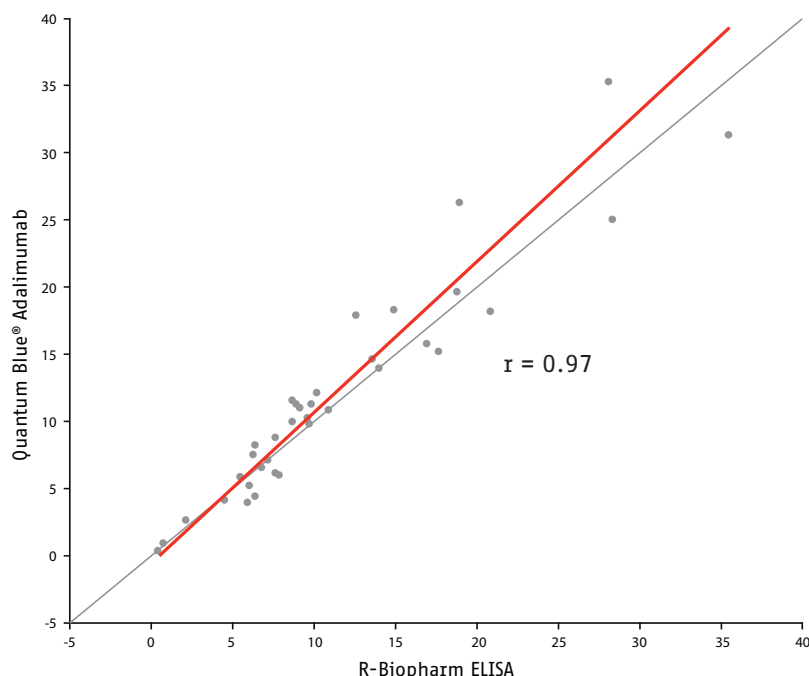

— Passing-Bablok fit

		RIDASCREEN ADM MONITORING (ELISA, µg/mL)		
		low	optimal	high
Quantum Blue® Adalimumab (µg/mL)	low	20.0%	5.0%	0.0%
	optimal	0.0%	27.5%	0.0%
	high	0.0%	12.5%	35.0%

total agreement 82.5%



BÜHLMANN Laboratories AG
Switzerland
info@buhlmannlabs.ch
www.buhlmannlabs.ch
Phone +41 61 487 12 12

BÜHLMANN Diagnostic Corp
USA / Canada
info@buhlmannlabs.com
www.buhlmannlabs.com
Phone 844 300 9799

Health Canada Licence: 101776
Quantum Blue® Adalimumab (LF-TLAD10/LF-TLAD25)

Quantum Blue® is a registered trademark of BÜHLMANN in many countries.

Ordering codes:
LF-TLAD10 10 cassettes
LF-TLAD25 25 cassettes

Related products:
LF-TLIF10 10 cassettes
LF-TLIF25 25 cassettes
LF-CAL25 25 cassettes
LF-CALE25 25 cassettes
LF-CHR25 25 cassettes